

Size and colour chart

Ankle High Sock – Adult Open or Closed Toe


SIZES

Size	Height from floor (a)	Top circumference (b)	Ankle circumference (c)	Foot length (d)
X SML	11cm	11 – 20.5cm	28.5 – 29.5cm	22 – 23cm
SML	11.5cm	20.5 – 22cm	29.5 – 30.5cm	23 – 24cm
MED	12cm	22 – 23.5cm	30.5 – 31.5cm	24 – 25cm
LGE	12.5cm	23.5 – 24.5cm	31.5 – 32.5cm	25 – 26cm
X LGE	13cm	24.5 – 25.5cm	32.5 – 33.5cm	26 – 27cm
XX LGE	13.5cm	25.5 – 27cm	33.5 – 34.5cm	27 – 28cm




MEASUREMENTS

Measuring Tips:




- All circumferences and lengths to be measured to the nearest millimetre.
- Tape should be snug NOT tight.
- Take (a) measurement to determine location for (b) measurement.

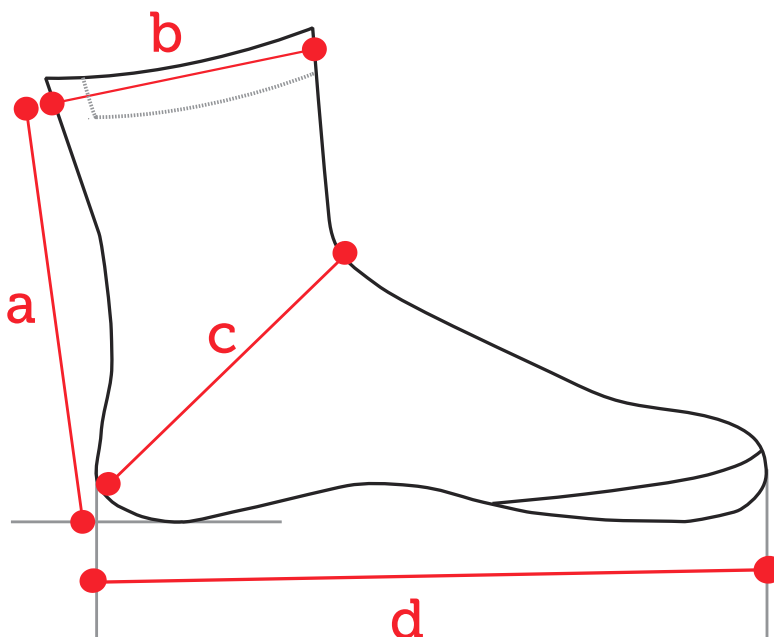

FABRIC OPTIONS

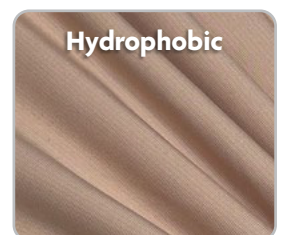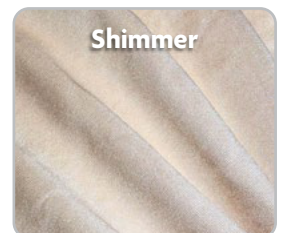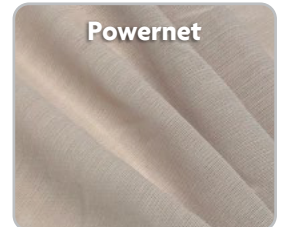

COLOURS

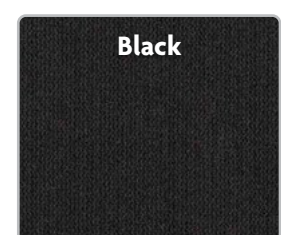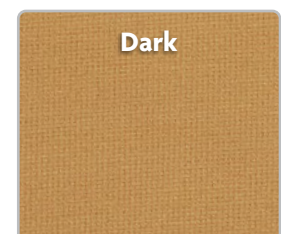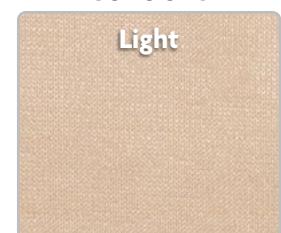

Please note garments *cannot* be altered or repaired by Second Skin.



Second Skin Pty Ltd
40 O'Malley Street
Osborne Park WA Australia 6017
Email: perth@secondskin.com.au
Phone: +61 8 9201 9455



Second Skin
10/32-34 Gordon House Road
Gospel Oak, LONDON NW5 1LP
Email: london@secondskin.com.au
Phone: +44 203 627 9937

secondskin.uk





Second Skin

CUSTOMISED COMPRESSION GARMENTS & SPLINTS

Fabric options

All Second Skin garments are made from a range of medical grade stretch fabrics. Several highly durable medical grade fabrics are incorporated in our garment designs, to address the individual needs of each client.

POWERNET FABRICS

Powernet is an open weave fabric which is prescribed for healed skin and for areas where additional support and compression is required. This fabric is particularly durable in work, school or outdoor activities. Powernet is often used in burns scar management and lymphoedema.

Our Powernet fabric now incorporates patented cooling properties, HeiQ Smart Temp, which provides outstanding thermal comfort to clients by keeping them cooler and dryer while wearing their garment.

Powernet fabric is available in light tan, dark tan and black.

SHIMMER FABRICS

Shimmer is a closed weave fabric which provides a softer interface with the skin. Shimmer fabric contours well and is comfortable and effective over joint surfaces and areas of skin fragility such as the elbow, ankle and knee. Shimmer fabric is used for early burns and flap management, external fixators, vascular conditions and certain types of lymphoedema.

Shimmer fabric is available in light tan, dark tan and black.

HYDROPHOBIC FABRICS

Hydrophobic fabric is used as both a lining and a base fabric in compression garments. This fabric wicks moisture away from the surface of the skin to evaporate, protecting the skin surface and ensuring optimal comfort.

Our hydrophobic fabric has an antimicrobial component embedded within the matrix of the fabric. This ensures our products are breathable and highly effective at removing bacteria from the fabric surface, to provide the wearer with all year round comfort. Due to this innovation, our garments are frequently prescribed to assist with the management of repeated skin breakdown and fragility.

Hydrophobic fabric is available in light tan, dark tan and black.

PERTH OFFICE

40 O'Malley Street
OSBORNE PARK
WA 6017

P: +61 8 9201 9455

E: perth@secondskin.com.au

UK OFFICE

Unit 10, Spectrum House
32-34 Gordon House Road
Gospel Oak, LONDON NW5 1LP

P: 0203 6279937

E: london@secondskin.com.au

 Follow us on LinkedIn